## Experiment



To identify washing soda or baking soda among given samples of chemicals.

# THEORY 💝



Washing soda (Na<sub>2</sub>CO<sub>3</sub>.10H<sub>2</sub>O) and baking soda (mainly NaHCO<sub>3</sub>) are white solids. Their aqueous solutions are alkaline and turn red litmus blue. Carbonates and hydrogenearbonates react with dilute acids and produce carbon dioxide gas which turns lime water milky.

$$CO_3^{2-}(aq) + H_2O(l) \rightarrow H_2CO_3(aq) + 2OH^-(aq)$$

$$HCO_3^-(aq) + H_2O(l) \longrightarrow H_2CO_3(aq) + OH^-(aq)$$

$$Na_2CO_3(s) + 2HCl(aq) \longrightarrow 2NaCl(aq) + H_2O(l) + CO_2(g)$$

$$NaHCO_3(s) + HCl(aq) \longrightarrow NaCl(aq) + H_2O(l) + CO_2(g)$$

$$Ca(OH)_2(aq) + CO_2(g) \longrightarrow CaCO_3(s) + H_2O(l)$$

$$Lime water \qquad Milky$$

On passing excess of  $CO_2$  through limewater, calcium hydrogenear bonate is formed. It is soluble in water and forms a colourless solution.

$$CaCO_3$$
 (s) +  $H_2O$  (l) +  $CO_2$  (g)  $\longrightarrow$   $Ca(HCO_3)_2$  (aq).

## MATERIALS REQUIRED



Samples: sodium carbonate (washing soda), sodium hydrogencarbonate (baking soda), ammonium chloride, sodium chloride etc., red litmus paper strips, freshly prepared lime water, dil. hydrochloric acid, five test tubes, a test tube stand, a boiling tube, a thistle funnel, a double bored cork, a delivery tube, and a glass rod.

#### **PROCEDURE**





- 1. Take about 1 g each of the four given samples separately in four boiling tubes and label them as A, B, C, and D.
- 2. Add about 5 mL distilled water in each boiling tube. Gently shake the contents of the tubes.
- 3. Put a drop of every salt's solution on separate red litmus paper strips, using a glass rod. (Wash the glass rod used for one sample before using it for other sample.) Note the change in colour of the litmus paper, if any in each case.
- 4. Add 1 mL of dil. hydrochloric acid in each test tube. Do you see any effervescence from any test tube? If ves, perform

the lime water test as

detailed below.

5. For performing lime water test, take the solution of test tube A in a boiling tube and set up the apparatus (delivery tube, thistle funnel etc.) as shown in Fig. 9.1.

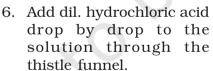

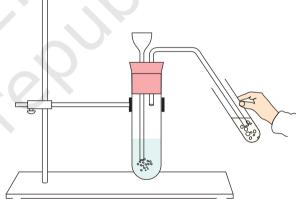

drop by drop to the Fig. 9.1: Carbon dioxide gas formed by the reaction of dil/hydrochloric acid on washing soda or baking soda is being passed through lime water

- 7. Pass the liberated gas evolved through the lime water in a test tube. Does the lime water turns milky? If yes, then it shows the presence of CO<sub>o</sub> gas.
- Continue passing the liberated gas through the lime water. Does it again become colourless? This reconfirms that the liberated gas is CO<sub>2</sub>.
- 9. Repeat the lime water test on all samples that give effervescence in step 4. Do not forget to wash the boiling tube when you change the sample in it for performing the lime water test.

#### **O**BSERVATIONS



Sl. No.	Sample	Colour	Solubility in water	Action on red litmus paper	Action of dil. HCl acid	Lime water
			(Soluble/ insoulble)	(Changes to blue or not)	(Effervescence observed or not?)	Turns milky ) or not?
1.	A					
2.	В					
3.	C					
4.	D					

#### RESULTS AND DISCUSSION



Infer from the observations about the identification of washing soda or baking soda out of the samples given for testing. Discuss about the litmus paper and lime water tests performed.

Sample in test tube \_\_\_ is washing soda/baking soda.

#### **P**RECAUTIONS



- Add dil. hydrochloric acid to the salt solution drop by drop. If the addition of dil. hydrochloric acid is not slow, a vigorous reaction may occur and the reaction mixture may come out of the reaction tube and pass into lime water.
- Handle hydrochloric acid and washing soda carefully. These should not touch your skin.
- Freshly prepared lime water should be used for performing lime water test

#### NOTE FOR THE TEACHER

- Students may be given three or four samples of salts, of which one of the salts is washing soda or baking soda. The remaining samples of salts should not be carbonates, hydrogenearbonates, sulphites or hydrogensulphites. These salts liberate either CO<sub>2</sub> or SO<sub>2</sub>. Sulphur dioxide also turn lime water milky.
- Preparation of lime water. Shake 5 g calcium oxide CaO, with 100 mL water. Allow it to stand for about 24 hours. Decant the supernatant liquid and use it for the tests. It is suggested to always use freshly prepared limewater.

### **Q**UESTIONS

- Explain why should dil. hydrochloric acid be added dropwise to the salt solution while performing lime water test?
- What will happen if crystalline washing soda is left open in air?
- CO<sub>2</sub> and SO<sub>2</sub> both turn lime water milky and their aqueous solutions turn blue litmus paper red. How can you then distinguish between these?
- Why should carbon dioxide be soluble in aqueous solution of potassium carbonate?